
Sirius TOMOGRAPH AND CORNEAL TOPOGRAPHER

Combines placido disk topography with Sheimpflug to- In addition to the clinical diagnosis of the anterior segmography of the anterior segment. Sirius provides infor- ment the most common uses are: refractive and catamation on pachymetry, elevation, curvature and dioptric ract surgery, an IOL calculation module is available. power of both corneal surfaces over a diameter of 12 Objective examinations provide an accurate measurment mm. All biometric measurements of the anterior cham- of pupil diameter in scotopic, mesopic and photopic conber are calculated using 25 sections from the cornea. ditions. When combined with the corneal map they can Measurement speed reduces the effect of eye move- be used for refractive surgery planning and follow up. ment producing a high quality accurate measurement.

FEATURES OF THE PHOENIX SOFTWARE

Sirius uses the Phoenix software platform allowing pa-On the basis of the pachymetry map and corneal altitient data to be saved for future review and analysis, shametric data, SIRIUS allows for intrastromal rings system red by all CSO devices. planning, which maybe an option for the correction of refractive defects and some forms of keratoconus.




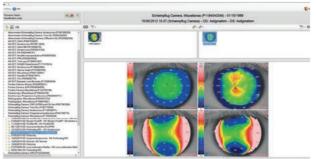




GLAUCOMA SCREENING

For glaucoma specialists Sirius enables the measurement of irido-corneal angles and pachymetry. These two values are useful in the diagnosis of the disease.

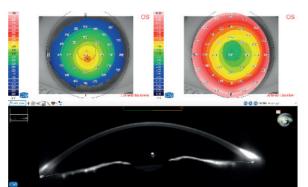


INTRASTROMAL RINGS

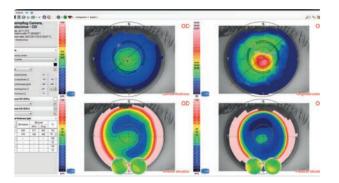


KERATOCONOUS SCREENING

Keratoconous screening provides the clinician with important information about the patinets cornea. Understanding this can help prevent complications associated with ectasia before corneal surgery is undertaken.

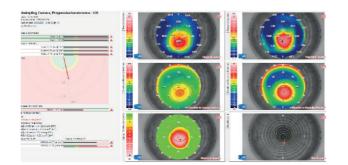


PUPILLOGRAPHY

sopic (4 lux), photopic (50 lux) conditions and in dynamic many lens manufacturers. mode. Knowing the center and the diameter of the pupil, is essential for many clinical procedures which seek to optimize vision quality.

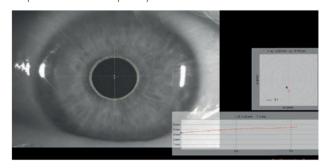


IOL CALCULATION MODULE (OPTIONAL)

This module is based on Ray-Tracing techniques, regardless of the state of the cornea (untreated or previously treated for refractive purposes), provides the calculation of the spherical and toric power of the intraocular lens.

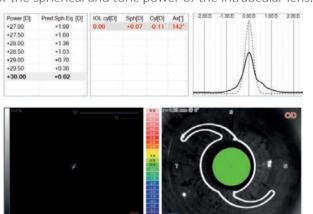


ADVANCED ANALYSIS OF THE TEAR FILM

Placido disk technology allows for the advanced analysis of the tear film, such as NI-BUT (Non Invasive Break-up Time).

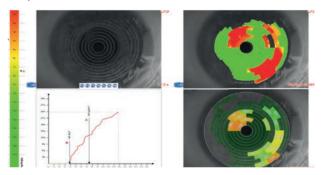


CONTACT LENSES APPLICATION MODULE

Sirius has built-in pupillography measurement software. A contact lens fitting module is available which simulates The measurement of the pupil in scotopic (0.04 lux), me- the fit of rigid lenses based on an internal database of

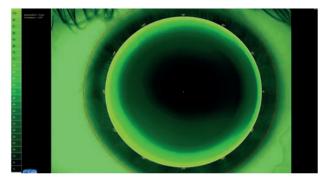

CORNEAL ABERROMETRY

Aberrometric analysis offers a complete overview of the corneal aberrations. It is possible to select the contribution of the anterior, posterior or total cornea for different pupil diameters. The OPD/WFE maps and the visual simulations (PSF, MTF, image convolution with optotype) can help the clinician in understanding or explaining the patient's visual problems.

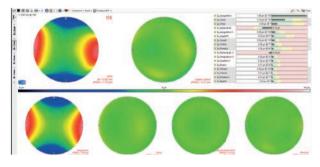

MEIBOGRAPHY

Meibomian glands can be viewed under infrared light once the image is captured, you can use the software to aid in the analysis of the condition of the glands



Sirius Tomograph and corneal topographer

TECHNICAL DATA


| Data transfer | Firewire |
|---|---|
| Power supply | external power source 24 VCC In: 100-240Vac - 50/60Hz - 0.9-05A - Out: 24Vdc - 40W |
| Power net cable | IEC C14 plug |
| Dimensions (HxWxD) | 515 x 315 x 255mm |
| Weight | 7 Kg |
| Chin rest movement | 70mm ± 1mm |
| Minimum height of the chin cup from the table | 24cm |
| Base movement (xyz) | 105 x 110 x 30mm |
| Working distance: | 74mm |
| LIGHT SOURCES | |
| Placido disk | Led @635nm |
| Scheimpflug | Led @475nm UV-free |
| Pupilligraphy | Led @875nm |
| TOPOGRAPHY | |
| Placido rings | 22 |
| Measured points | 21632 (front surface) 16000 (rear surface) |
| Topographic covering | 12mm |
| Dioptric measurement range | 1D to 100D |
| Measurement accuracy | Class A according to UNI EN ISO 19980-2012 |
| Compatibility with standard | DICOM v3 (IHE integration profile EYECARE Workflow) |

MINIMUM SYSTEM REQUIREMENT

PC: 4 GB RAM - Video Card 1 GB RAM (not shared) resolution 1024 x 768 pixels - Firewire port. Operating system: Windows XP, Windows 7 and Windows 10 (32/64 bit).

*The specifications and the images are not contractually binding and can be modified without notice. Windows® is a Microsoft Corporation trade mark.

CO130 | Rev. 00 del 01/2018

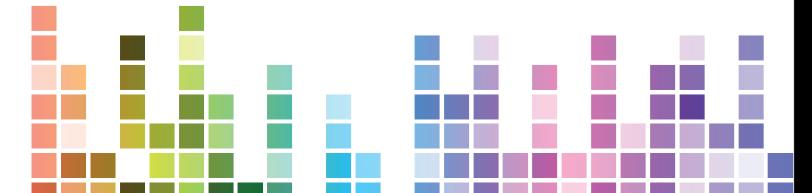





Via degli Stagnacci 12/E 50018 - Scandicci - FI - Italy tel +39 055 72219 | fax +39 055 7215557 email. cso@csoitalia.it | web. www.csoitalia.it



